



LIFT MB750 MICRO MINIMUM DIMENSIONS

Ultra slim pass-through lift

With the Herkules wheelchair lifts MB750 "micro" and "duo" you have two variants to get access to buildings with limited available space. The MB750 "micro" is designed as an ultra slim pass-through lift. The external dimensions can be reduced to only 94×136 cm, thanks to the proven Herkules spindle drive technology from the wheelchair lift MB1000 / MB1300 and the flat scissor- system. The Hercules MB750 "duo" makes a pass-through and diagonal entrance possible with only 125×150 cm footprint. The maximum adjustable lifting height of 75 cm covers up to four steps inside and outside the house. The MB750 lifts up heavy electric wheelchairs easily, soft and noiseless with its high loading capacity of 350kg. No matter if it is cold, hot or rainy – the weatherproof Mb750 is a lift solution you can rely on.



LIFT MB750 "MICRO"

LIFTING HEIGHT: 750 MM

LOAD: 350 KG

DIMENSIONS: 940 X 1.360 MM


CLEARANCE: 750 MM

INDOOR AND OUTDOOR USE

PIT DEPTH: ONLY 55MM

HERKULES REHA GMBH FALDERBAUMSTR. 36 34123 KASSEL

www.herkules-reha.de info@herkules-reha.de Weitergabe sowie Vervielfältigung dieses Dokuments, Verwertung und Mitteilung seines Inhalts sind verboten, soweit nicht ausdrücklich gestattet. Zuwiderhandlung verpflichtet zu Schadenersatz. Alle Rechte für den Fall der Patent-, Gebrauchsmuster- oder Geschmacksmustereintragung vorbehalten. Herkules Reha GmbH, Kassel, Deutschland

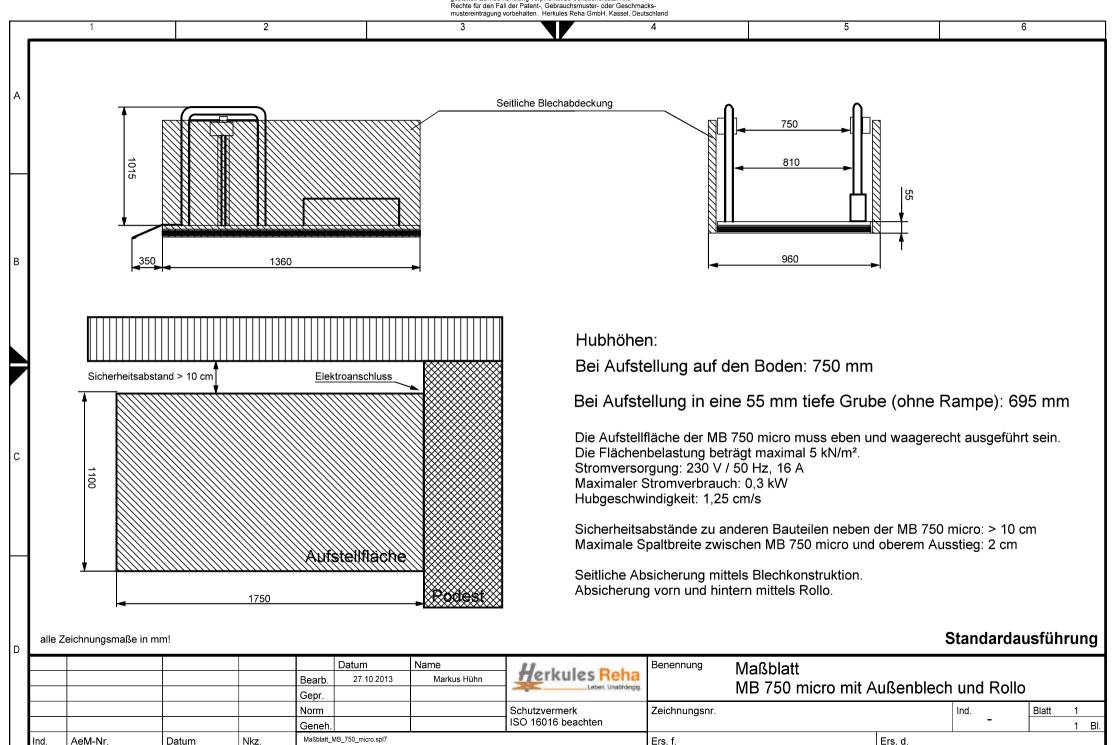

Weitergabe sowie Vervielfältigung dieses Dokuments, Verwertung und Mitteilung seines Inhalts sind verboten, soweit nicht ausdrücklich gestattet. Zuwiderhandlung verpflichtet zu Schadenersatz. Alle Rechte für den Fall der Patent, Gebrauchsmuster- oder Geschmacksmustereintgung vorbehalten. Herkules Reha GmbH kassel Deutschland

